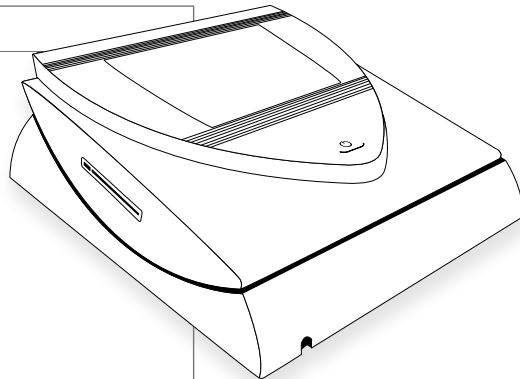
DEUTSCH
ENGLISH
FRANÇAIS
MAGYAR

BEMER-SET CLASSIC/PRO

:: User Manual

FOR YOU :: FOR LIFE :: FOR ENERGY

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Congratulations on the purchase of your new BEMER set. These high quality products are manufactured by us with the greatest of care and meet our high quality standards as well as the European and international production standards of medical products. Not only this high level of quality but also a high level of user-friendly functionality account for more than 10 years of well above average customer satisfaction.

MANUFACTURER:

BEMER INT. AG

Schliessa 12 // 9495 Triesen
Liechtenstein

WWW.BEMERGROUP.COM

NOTIFIED BODY:

LGA INTERCERT GMBH

Certification body of medical devices
Tillystraße 2 // 90431 Nürnberg
Germany

FOR YOU :: FOR LIFE :: FOR ENERGY

Contents

Scope of delivery	76	B.BOX Professional – display	110
Foreword	78	Applications – using the B.BOX Professional	112
Intended use	79	Set up – B.BOX Professional	120
Safety instructions	80	Updates	122
User information	84	B.SCAN	124
Examples of application	94	Rechargeable battery – B.BOX	126
Connections – B.BOX	98	Wall mount	128
Start-up procedure – B.BOX	99	Care / maintenance	130
Connection – application module	100	Proper disposal	131
General information	101	Troubleshooting	132
B.BOX Classic – display	102	Technical data	134
Applications – using the B.BOX Classic	103	Manufacturer's declaration of conformity	138
Set up – B.BOX Classic	108	Guarantee	144

Scope of delivery Classic Set* [Art. No. 410100]

Art. No. 420100

B.BOX Classic (control device with multifunctional display)

APPLICATORS

B.BODY Classic (application module for full body therapy)		[Art. No. 430100]
B.GRIP (holder module for B.SPOT/B.LIGHT)		[Art. No. 431000]
B.SPOT (application module for targeted therapy)		[Art. No. 431100]
B.SCAN (scanning module for functional checks)		[Art. No. 450100]

INCLUDED ACCESSORIES

Power pack B.BOX		[Art. No. 440100]
Attachment belt (B.GRIP)		[Art. No. 450700]
Attachment belt (B.BODY)		[Art. No. 450800]

* Please check the scope of delivery for completeness.
Please note which set you have purchased.

Scope of delivery Pro-Set* [Art. No. 410200]

Art. No. 420200

B.BOX Professional
(control device with multitouch display)

APPLICATORS

B.BODY Professional (application module for full body therapy)		[Art. No. 430200]
B.GRIP (holder module for B.SPOT and B.LIGHT)		[Art. No. 431000]
B.SPOT (application module for targeted therapy)		[Art. No. 431100]
B.LIGHT (application module for light therapy)		[Art. No. 431200]
B.PAD (flexible application module for small surface treatment)		[Art. No. 430300]
B.SCAN (scanning module for functional checks)		[Art. No. 450100]

INCLUDED ACCESSORIES

Wall mount	[Art. No. 450200]	Foot protection	[Art. No. 450500]
Protective goggles	[Art. No. 450600]	Power pack B.BOX	[Art. No. 440100]
Attachment belt (B.GRIP)	[Art. No. 450700]	Vehicle connection cable	[Art. No. 440300]
Attachment belt (B.BODY)	[Art. No. 450800]		

Foreword

Symbols found in the operating instructions

Warning / Caution – failure to observe can result in possible injuries or property damages.

Eye protection – failure to observe can result in possible injuries to the eyes.

Proper disposal

Tip

Manufacturer

CE conformity marking according to Medical Devices Directive 93/42/EWG

The identification number of the Notified Body LGA Intercert GmbH

The BEMER therapy is designed to support the body's fundamental physiological processes, thereby promoting the improvement of natural self-regulatory mechanisms and offers the following scientifically proven effects: positive effects on microcirculation, the immune system, the protein biosynthesis and the formation of free radical scavengers that have been naturally produced by the body.

Therefore, the BEMER therapy already offers a wide range of application possibilities with very good supportive effects in the areas of therapy, regeneration, rehabilitation and prevention.

BEMER can be used as a supplementary application wherever disturbances of the microcirculation occur that have been caused by the progress of diseases or by accompanying illnesses, or that have led to deficiency symptoms.

Intended Use

By means of a weak, pulsating and special multidimensional signal configuration of generated electromagnetic fields, BEMER products are suitable for stimulating vasomotion (rhythmic contractions) of the smaller and smallest blood vessels in living organisms. Through certain biorhythmic modulations, somewhat larger blood vessels upstream and downstream are also synergistically influenced. This results in more favorable microcirculatory blood flow for the treated organisms. The results are: improved functionality of natural regulatory mechanisms, strengthening of the immune system, improved production of the body's own radical scavengers and an influence on the protein synthesis. The BEMER technology can be used for the prevention and complementary treatment of all diseases, leading to an improvement in microcirculation and resulting in a definite therapeutic effect in living organisms.

The application module B.LIGHT complements the user benefits of the BEMER therapy by means of photobiological light mechanisms. Through its light effect, the B.LIGHT promotes the regeneration processes of the skin. Extremely bright light-emitting diodes (LED) emit a pulsating, monochromatic red light with wavelength of 660 nm, making the B.LIGHT an ideal supporting application for skin diseases and also for supplementary cosmetic application.

BEMER products may not be used for any purpose other than that described in this chapter. Use of the device beyond that which has been specified is considered contrary to the intended purpose. The manufacturer is not responsible for any damages resulting from improper use. The user is solely responsible for any risks.

Safety instructions

Device symbols	
CE conformity marking according to Medical Devices Directive 93/42/EWG	
The identification number of the Notified Body LGA Intercert GmbH	1275
Proper disposal	
Manufacturer	
Please observe operating manual	
Protection class Device is of the protection class II	SK II 
Application part type Application part is the type BF	BF 
Serial number	

The BEMER products are manufactured according to the latest state of technology and are completely safe to operate. Nevertheless hazards could emanate from this device. This is particularly the case when insufficiently trained personnel are operating the device or should the device be inappropriately used and not according to its intended use.

- >> Complaints, symptoms or diseases must be clarified by a doctor / therapist before application with the BEMER technology begins, should this not already be the case.
- >> Prescribed medications may only be altered after consulting with the treating physician.
- >> Unexplainable reactions of the body, which occur during the BEMER therapy, are to be discussed with the physician in charge, for his / her appropriate treatment if necessary, in the event of an emerging disease or illness.
- >> Application with BEMER systems is not allowed for persons with organ transplantations, allogenic stem cell transplantations or bone marrow transplantations.

-
- >> The therapeutic use of BEMER technology is not allowed for persons who have foreign organs or have had cell transplants and therefore have been medically prescribed with immunosuppressants.
 - >> Persons that are affected by infectious disease, severe cardiac arrhythmia, severe psychosis and seizures that are not compensated (e.g. epilepsy) must consult with their physician in charge before beginning application with the BEMER technology.
 - >> Persons with electronic implants may only use BEMER technology after confirming that the limits set by the implant manufacturer are not exceeded. Interactions or interferences with electronic implantations such as e.g cardiac pacemakers, ICDs, brain pacemakers, insulin pumps, etc. can not be generally excluded.
 - >> The BEMER therapy is an important support therapy for physicians / therapists and may not replace the treatment of doctors / therapists.
 - >> Before you connect the BEMER devices to the power supply, please make sure that the power supply is identical with that of the B.BOX power pack supply.
 - >> Power cable and power plug must be in their original, undamaged condition. Immediately replace any damaged parts!
 - >> Make sure to unplug the power plug when cleaning!
 - >> Do not misuse BEMER devices! The devices may not be used for any purposes other than those mentioned here.
 - >> Always keep the BEMER devices in proper working order!
A neglected device increases the risks of accidents.
 - >> Only IECxxxxx certified components may be connected to the maintenance / communication interface.

- >> Please check the BEMER devices on a regular basis for damage. Immediately replace any damaged parts!
- >> Please observe the user information.
- >> BEMER devices may only be used on a stable, non yielding surface (e. g. floors, sofa, bed, table).
- >> BEMER devices are not toys.
Keep children away from devices.

Special safety instructions for the B.LIGHT

- >> It is recommended that all skin diseases caused by bacteria or fungi as well as light allergies and phototoxic reactions be treated by a physician.
- >> Light emitted by the LEDs may not be pointed directly into the eyes.
Use the provided protective goggles!
- >> Make sure that no third person looks directly into the light or is illuminated by the LEDs.
- >> Do not use any sharp edged or pointed objects to clean the B.LIGHT.

MAGYAR

FRANÇAIS

ENGLISH

DEUTSCH

User information

These operating instructions contain important information. Keep them close to your device in order to inform yourself about safety instructions as well as proper handling.

Please read the operating instructions very carefully before using the device. In this way you can be assured that you benefit from all advantages that this system has to offer and protect yourself and others from harm.

The BEMER devices are designated only for the described, intended use according to the operating instructions. In the case of foreign translations: the original German version is binding.

Application principles

1. GENERAL FULL BODY THERAPY

The standard general full body therapy is carried out with the application module B.BODY* using a maximal intensity of 35 μT that corresponds to the basic plan (BP) and is used for the general improvement of the overall physical well being of the body. It improves the requirements necessary for varied metabolic processes and self-regulating mechanisms. This forms the basis for its application in prevention and therapy.

B.BODY is the standard application module used for full body therapy and has the same positive effect regardless of which side of the mat you choose to use. The body position on the B.BODY* should be in the middle of the mat. It is unimportant if the head, legs or both of taller users extend beyond the mat ends. You will recognize the head area of the B.BODY from the connection cable.

For use with the Sleep Program, the B.BODY can be installed in the bed with the help of the attachment belt.

People with limited physical movements are “covered” by the B.BODY*, which they can use in either a sitting or lying position.

As countless applications in therapeutic practices have shown, people react very differently to the individual intensity levels, in which numerous correlations between metabolic status, intensity and effect are recognized. When using the BEMER Signal configuration, different reactions of circulation parameters at various intensity levels have been observed.

2. TARGETED APPLICATION

The additional application modules B.SPOT and B.PAD are intended for targeted and intensive application of individual body parts, working with a maximum intensity of 100 μ T. Basically, the targeted application should follow the general application, since important fundamental processes are set into motion with the targeted application can now optimally build upon. Multiple applications are possible in one day.

B.SPOT is an intensively active application module for use in localized and targeted therapy. Because of its small size it's perfect for travel. The B.SPOT is attached to the universal B.GRIP holder and held on the body area that is to be treated. Alternatively, the holder and application module can also be attached to the part of the body that is in need of treatment with the provided attachment belt.

B.PAD is an intensively active application module for targeted, small surface application on individual body parts. The B.PAD is especially attractive because of its flexibility and can be attached with the help of a Velcro fastener strap to nearly every conceivable part of the body. When using on disabled or mobility impaired people, the B.PAD can provide a particularly valuable service.

B.LIGHT is an application module used for complementary therapy of skin diseases, wound healing disorders and for cosmetic application. The light-emitting diodes create a monochromatic and highly effective photobiological red light. This stimulates the skin and connective tissue cells to produce additional energy-rich compounds. The effect of the B.LIGHT is comparable to the effect of soft laser treatments.

Therefore, high intensity levels are usually recommended. It is also possible to select lower intensities when using on very sensitive skin. The complementary and the subsequent full body application is recommended for, among other things, acne, eczema, infantile eczema, neurodermatitis, psoriasis, skin injuries or warts.

3. PROGRAMS AND INTENSITY LEVELS

The control devices B.BOX Classic and B.BOX Professional offer 10 different intensity levels and three predefined programs. The intensity levels are applied during the generalized full body application according to the basic plan while the programs P1 – P3 are carried out in a step-by-step “deep penetration” within the scope of the targeted application. The existing, predefined programs cover the wide range of application requirements needed for daily use.

All the same, it is possible for you to create and save your own programs (only with the B.BOX Professional). Application programs prepared by you, your doctor or therapist can be stored on a memory card and performed on every B.BOX.

Frequency of application

Whenever possible, the BEMER therapy should be carried out on a daily basis. Intermittent applications mean a lessened effect, resulting in longer application durations. Overdosing or habituation effects of the BEMER system are not known.

There is no hard and fast rule for the distribution of applications throughout the day. It depends on the therapy goal and the situation of the individual user.

- >> For prevention and for maintaining health and the body's performance, one to two applications during the day, as well as carrying out the Sleep Program, are perfectly adequate.
- >> Occasional frequent or consecutive treatments using the application modules may be necessary in the event of acute injuries.

- >> At home application can be customized according to the individual's situation and needs. If possible, applications should be carried out once in the morning and again in the evening.
- >> At the doctor's office or therapeutic practices, it is usually possible to have only one application per day. This generally consists of one application with the B.BODY and immediately following, another using an application module that is especially designed for the purpose of the localized application (double application).
- >> Sometimes however, certain time constraints allow only one application a day. In this case, treatments would be similar to those at a doctor's office or therapeutic practice, where double applications are recommended.

Control device programs

The BEMER therapy is based on the nature of illness and on reduced performance and is not a classic indication-related treatment. Therefore it is not necessary to include an enumeration of medical conditions.

PROGRAM	DURATION	PENETRATION DEPTH IN THE ORGANISM
P1	8 min.	Low intensity for superficial / minor symptoms
P2	16 min.	Middle intensity for somewhat deeper / moderate symptoms
P3	20 min.	Strong intensity for deeper lying / severe symptoms

Practical tips for application

The effects of the multidimensional signal configuration are scientifically proven and essentially independent of external influences. Nevertheless, the effects of the application can be optimized by taking into account these few tips:

- >> Make sure you find a comfortable position.
- >> Avoid wearing tight or constricting clothing.

- >> Choose a room with a comfortable temperature for your therapy.
- >> Avoid stress or distraction.
- >> Abstain from coffee or tea before application.
- >> Refrain from smoking cigars, pipes or cigarettes before application.
- >> Drink a glass of sparkling free mineral water before and after application.

Basic plan

The basic plan is targeted at the optimal support of regulatory mechanisms with the corresponding positive physiological effects. It is the basic treatment for prophylaxis and for all pathological or any limiting conditions of a person's well-being, as well as for any unclear metabolic status.

The basic plan application (BP application) using the B.BODY* should generally be carried out twice daily and once while asleep (Sleep Program). Combinations of the B.BODY* with other levels, programs or using other application modules (double application) are possible.

The BP application cycle lasts 4 weeks. The application begins with level 3 and increases one level each week up to level 6. In the fifth week, a new cycle begins again using level 3.

Both the basic plan and all other application suggestions are based on many years of experience. This experience has also proven that regular applications are more important than the optimal choice of intensity. Successful therapy results are achieved quickest through regular application, however the duration of the applications depends on how difficult and how long the body's own regulatory systems have been impaired.

BASIC PLAN (Application using the B.BODY)

	WEEK 1	WEEK 2	WEEK 3	WEEK 4	WEEK 5
MORNING	Level 3 <i>plus</i>	Level 4 <i>plus</i>	Level 5 <i>plus</i>	Level 6 <i>plus</i>	Level 3 <i>plus</i>
EVENING	Level 3	Level 4	Level 5	Level 6	Level 3
NIGHT	Sleep Program	Sleep Program	Sleep Program	Sleep Program	Sleep Program

Sleep Program

The latest research results have proven that while sleeping, a redistribution of blood in the organism takes place and the immunological processes during this time have far more importance than previously assumed. Therefore, in addition to daytime applications, a further special signal configuration has been developed, which takes into account the physiological conditions that take place during sleep.

These regeneration processes during sleep are considered of major importance, particularly for rehab patients and the elderly, or for people susceptible to infections.

Before going to bed, the Sleep Program is activated on the control device and then automatically shuts itself off.

Advantages of the new Sleep Program

- >>A general increase in the application effect.
- >>Reduced susceptibility to infection.
- >>Greater stimulation of regenerative and restitutive (restorative) processes.
- >>A higher excretion of urinary waste products and acidic substances.

Medical hotline

Do you have any questions or need something clarified?
Then call our medical hotline and take advantage of our
professional advice.

Further information can be found at:

WWW.BEMERGROUP.COM

Examples of application

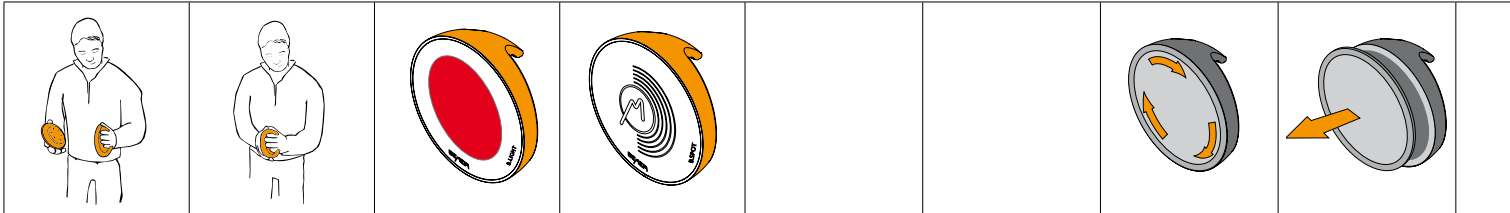
B.SCAN

The B.BOX has an integrated test device for checking the performance of your application module and can also be used for giving effective demonstrations.

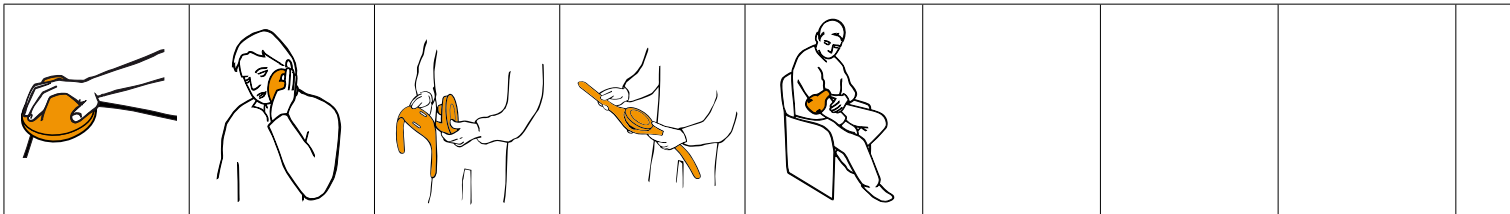
B.BODY Classic / B.BODY Professional
Application module for full body therapy.

B.GRIP

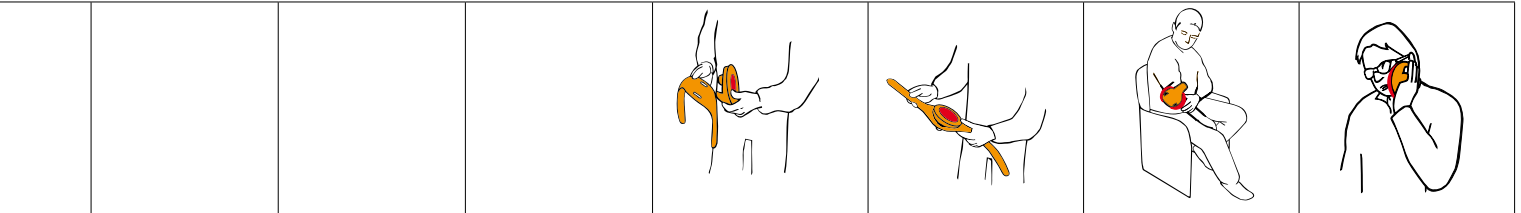
Multi-functional grip for B.SPOT and B.LIGHT.

B.SPOT

Flexible, portable application module for targeted application.

Examples of application

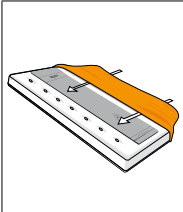
B.LIGHT

The BEMER Signal is an innovative, supportive application module with complex effects, for treatment using photobiologically effective light.

B.PAD

Handy, flexible and universally mobile. This application module is made for small surface treatment of specifically targeted areas of the body.

Fixing belt (B.BODY)

For attaching the B.BODY onto all standard mattresses and supplying additional hygiene.

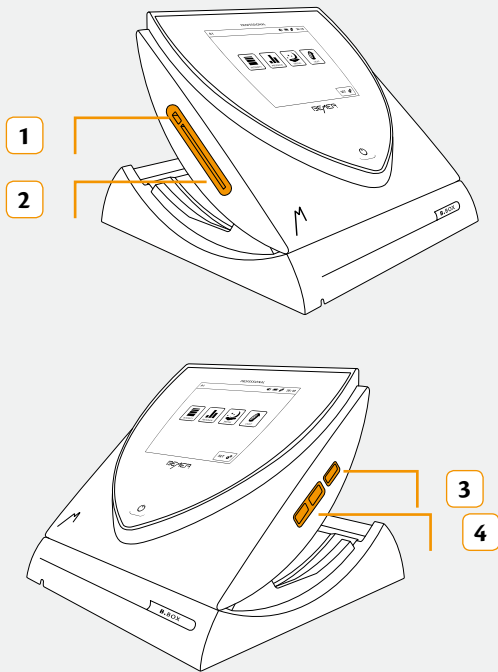
Made from extremely stretchable and light weight material.

Wall mount

A wall mount for easy fixation of the B.BOX –
for example, on the wall next to a double bed.

Connections – B.BOX

Maintenance interface¹

1

Memory card slot

2

Communication interface^{1 2}

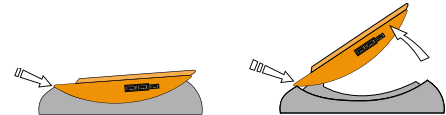
3

Connections for the application module (A1 and A2)

4

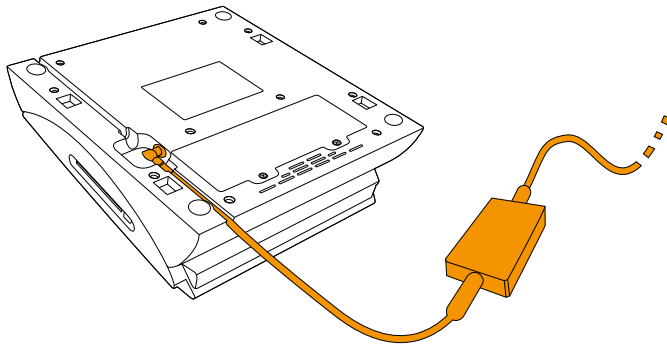
Basic position, recommended when using the wall mount

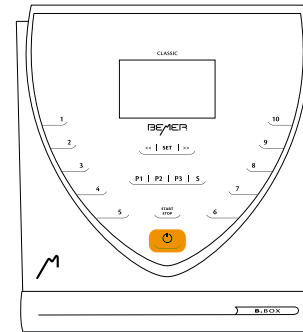
Position 1 and position 2, recommended for table application with different angles.

- ¹ Only IECxxxxx certified components may be connected to the maintenance / communication interface. Suitable for software updates and approved BEMER accessories.
- ² Applies only for B.BOX Professional.

Start-up procedure – B.BOX

Turn the B.BOX over and connect it to the adapter. Place the device on a stable surface and plug the adapter into the power supply.

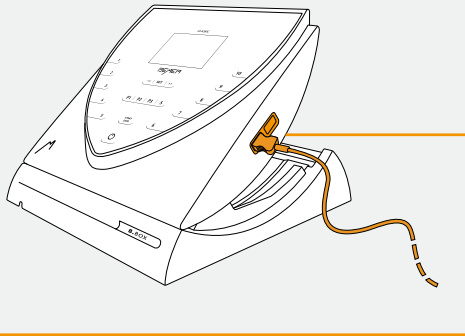
A short press on the power button switches the device on, and a longer press will switch it off.

After 2 minutes of inactivity,
the device automatically switches off.

Connection – application module

The magnetic plug prevents incorrect connection of the application module.

Connection of application modules

The two ports, A1 and A2, are available for connecting the application module. These are not specifically assigned, which means that they can connect any application module with any port.

General information

There is an important difference between the general and targeted forms of application when using BEMER therapy. The general full body application for day is carried out using the B.BODY according to the basic plan with the intensities 1-10 (short application). Application during the sleep phase is also carried out using the B.BODY and with the corresponding special program.

The targeted applications are carried out using the application modules B.SPOT and B.PAD and use, depending on the required penetration depth, the predefined programs P1 through P3.

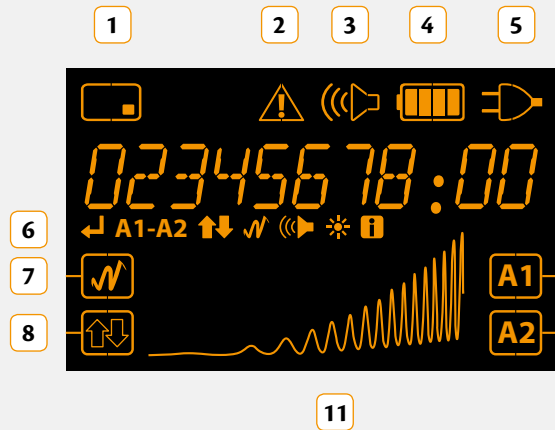
General tips

>> Please note that several application modules can be connected to the B.BOX.

Simultaneous operation when using the B.BOX Classic is only possible with two B.BODY application modules (using the same duration and intensity level).

- >> Short applications always last 8 minutes, regardless of the chosen intensity.
- >> Programs have different durations and use varying intensities.
- >> Ongoing applications can be canceled at any time by pressing the power button [START / STOP].
- >> Canceled applications can not be continued; they must be restarted.
- >> The *plus* Signal is switched on by default and can be switched off if needed.
- >> The Sleep Program does not support the *plus* Signal.
- >> It is not recommended to use the Sleep Program with the rechargeable battery.

B.BOX Classic – display

Memory card (if inserted)	1	<i>plus</i> Signal (if active)	7
Troubleshooting (if available)	2	Signal intensity (if active HIGH / LOW, only with B.PAD)	8
Volume (in current setting)	3	Info line	9
Rechargeable battery (state of charge)	4	Applicator (lit when connected, blinks when selected)	10
Power supply operation (if plugged in)	5	"Signal" (with ongoing application)	11
Menu navigation "Settings"	6		

Applications – B.BOX Classic

Short application (intensity levels 1-10)

- >> Chose the desired intensity by pressing the numbers 1-10.
- >> Designate the desired application module [A1 / A2] by repeatedly pressing [SET].

- >> Activate / deactivate the *plus* Signal by pressing the left arrow key.
- >> Start the short application with [START / STOP].

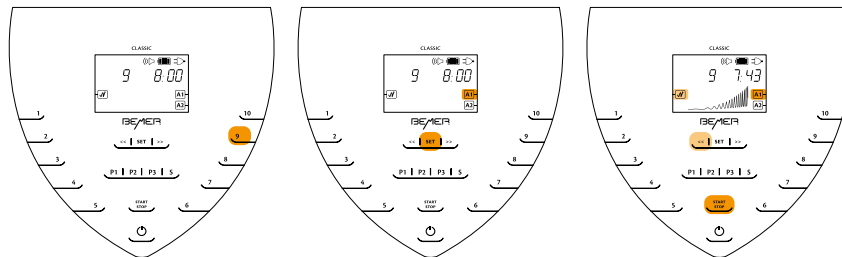
After selecting the desired application, press the left arrow key to switch the *plus* Signal on / off.

Simultaneous operation is only possible with 2 B.BODY application modules connected (both using the same duration and intensity).

Example – short application with intensity 9 (2 B.BODY application modules connected, port A1 selected)

The selected application module blinks both during the selection and during application use.

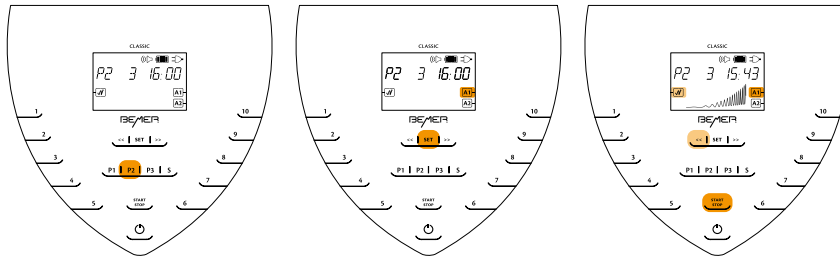
Program selection

- >> Select the desired program by using the keys [P1], [P2] or [P3].
- >> Designate the desired application module [A1 / A2] by repeatedly pressing the key [SET].

- >> Activate / deactivate the *plus* Signal by pressing the left arrow key.
- >> Start the application with the key [START / STOP].

Example – Program P2 (2 B.BODY application modules connected, port A1 selected)

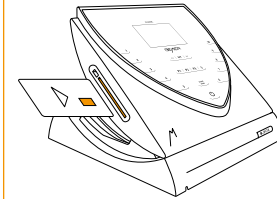
Application using the memory card

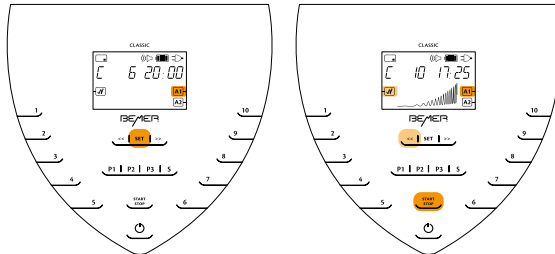
- >>Slide the memory card into the card slot.
- >>Designate the desired application module [A1 / A2] by repeatedly pressing the key [SET].

- >>Activate / deactivate the *plus* Signal by pressing the left arrow key.
- >>Start the program with [START / STOP].

The card must be removed again for short or program applications. 

Example – application using memory card (2 B.BODY application modules connected, port A1 selected)

Special application – Sleep Program

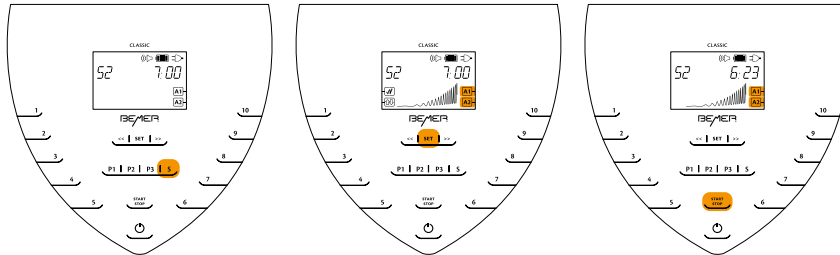
- >> Activate the Sleep Program by pressing the [S] key.
7:00 appears on the display (= 7 hours).
- >> Select the desired duration (6, 7 or 8 hours).

- >> Designate the desired application module [A1 / A2] by repeatedly pressing [SET].
- >> Start the program with [START / STOP].

The display screen fades after 90 seconds, so as not to disrupt your sleep. It can be reactivated at any time by pressing any key once.

Example – Sleep Program (2 B.BODY application modules connected, both simultaneously in operation)

Light therapy

- >> Designate the application module [A1 / A2], by repeatedly pressing the [SET] key.
- >> Select the desired light intensity by selecting a number [1] through [10].
- >> Start the program with [START / STOP].

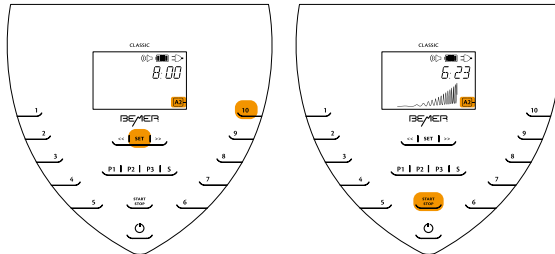
Usually high intensity levels are recommended. However for those with very sensitive skin, lower intensity levels can also be selected.

Do not let the LED light shine directly in the eyes. Use the provided protective goggles!

Example – using light therapy (B.LIGHT connected on port A2)

Default settings – application parameters / basic settings– B.BOX Classic

The parameters [PLUS] and [INTENSITY] can be adjusted to become the default setting in the settings menu. To access this menu, press and hold down the [SET] key. You can select the parameters [PLUS] or [INTENSITY] by using the arrow keys. To exit the menu, select with the arrow keys [EXIT] and confirm with the [SET] key.

Setting the *plus* Signal as default

- >> In the settings menu, select the menu item [PLUS] with the arrow keys and confirm with [SET].
- >> By using the arrow keys, the *plus* Signal can be switched on / off.
- >> Once the selection has been made, confirm with [SET]. Now the plus Signal has been switched on or off.

Default setting signal intensity (only with the B.PAD)

While traveling, the B.PAD can be used as a substitute for the B.BODY by reducing the intensity to [LOW].

- >> Select the menu item [INT] in the settings menu using the arrow keys and confirm with [SET].
- >> By using the arrow keys you can choose between [HIGH] (high intensity) and [LOW] (low intensity).
- >> Once the selection has been made, confirm with [SET].

Setting the signal volume

- >> Press and hold down the [SET] key to activate the settings menu.
- >> With the arrow keys, select [VOLUME] and confirm with [SET].
- >> Then choose the desired signal volume by using the arrow keys.

The selection will always be acknowledged with a signal in the chosen volume.
- >> Confirm selection with [SET].

Setting the display brightness

- >> Press and hold down the [SET] key to activate the settings menu.
- >> Using the arrow keys select [Display] and confirm with [SET].

>> Set the desired brightness by using the arrow keys.

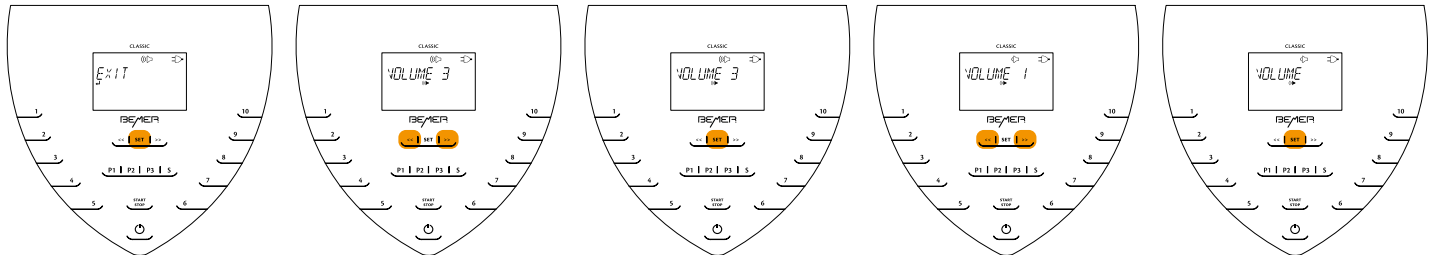
The display brightness corresponds to the current selection.

>> Confirm selection with [SET].

To access the system information

- >> Press and hold down the [SET] key to activate the settings menu.
- >> Select the menu item [INFO] by using the arrow keys and confirm with [SET].
- >> Now continue to your selection using the arrow keys.

Example – settings VOLUME

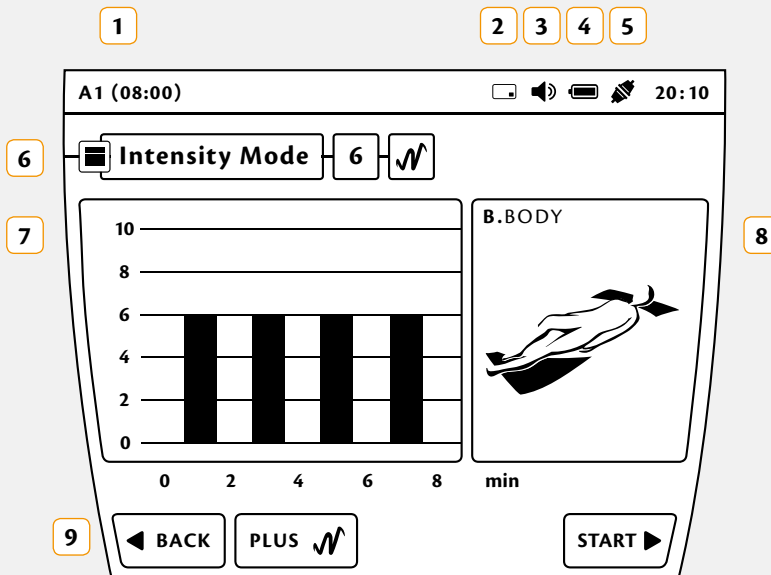
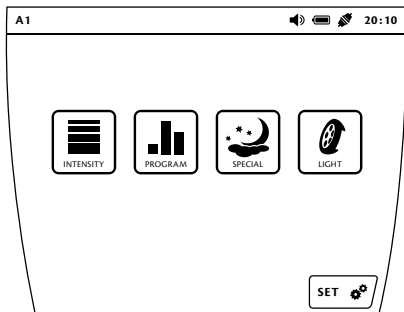
B.BOX Professional – display

Wait several seconds after turning on the device unit until it is ready for operation (see display below).

From the main window you can access the individual menu items by pressing the various buttons.

Applicators (display the remaining time of the ongoing application)

1

Memory card (if inserted)

2

Volume
(in current setting)

3

Rechargeable battery (state of charge)

4

Power supply operation
(if connected)

5

Displays the active modus
(the intensity and status
of the *plus* Signals)

6

Application window
(displays the application progress
incl. level and period of time)

7

Applicator window (displays the connected application module on the chosen port)

8

User and navigation lines (contain important operating elements, e. g. *plus* Button or [START / STOP])

9

Applications – using the B.BOX Professional

2-in-1 functionality. The absolute highlight of the B.BOX Professional is the possibility to simultaneously control and operate 2 of the same or even different application modules independently.

By pressing the buttons [A1 / A2] you can switch between the connected application modules, whereby it is possible, for example, to use the B.PAD and the B.LIGHT at the same time.

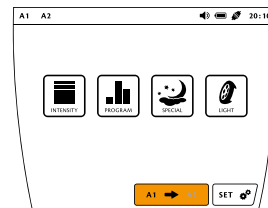
A further highlight of the B.BOX Professional is the option to operate 2 B.BODY Professional modules with different wake up times in the Sleep Program.

Example

A B.PAD is connected on port [A1] and a B.LIGHT on port [A2] and you want to combine a short application with light therapy.

- >> Press the button [A1 / A2] and proceed as described on page 41.
- >> Now press the button [A1 / A2] and proceed as described on page 45.

The remaining duration of both applications will be displayed in the upper information bar.

The [A1 / A2] button is only available when 2 application modules are connected.

Short application (intensity levels 1-10)

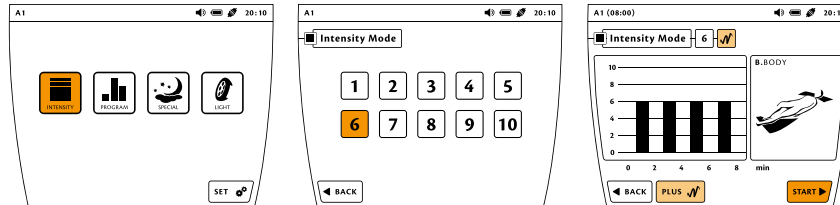
>> Press the button [INTENSITY].

>> Select the desired intensity by pressing the respective number (1-10).

>> Activate or deactivate the *plus* Signal by pressing [PLUS].

>> After making the selection, start the short application by pressing [START].

Example – short application with intensity 6 (using the application B.BODY via port A1)

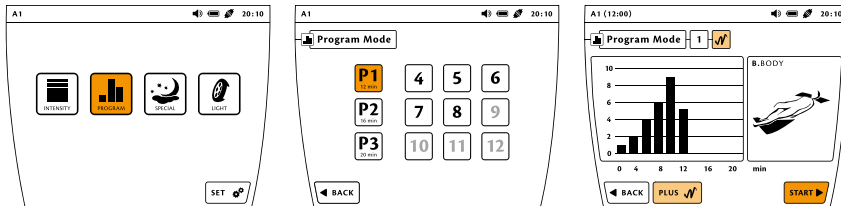
Program selection

- >> Press the button [PROGRAM]; if a memory card is inserted, press the button [PROGRAM] again.
- >> Select the desired program by pressing the button of one of the predefined programs P1 – P3 and as needed, the available programs 4 – 12.

- >> Activate or deactivate the *plus* Signal by pressing the button [PLUS].
- >> After making the selection, start the short application with [START].

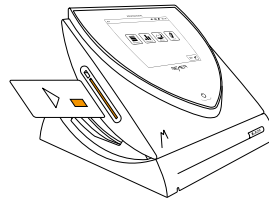
Example – Program P1 (with the application module B.BODY via port A1)

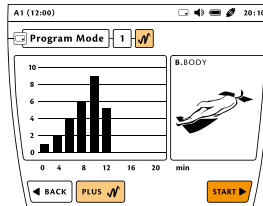
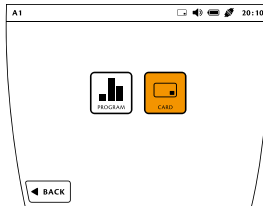
Application using the memory card

- >> Slide the memory card into the card slot.
- >> Press the button [PROGRAM].
- >> Press the button [CARD].
- >> Start the program with [START].

The selection [CARD] is available only when a memory card is inserted.

Example – application using a memory card (using the application module B.BODY via port A1)

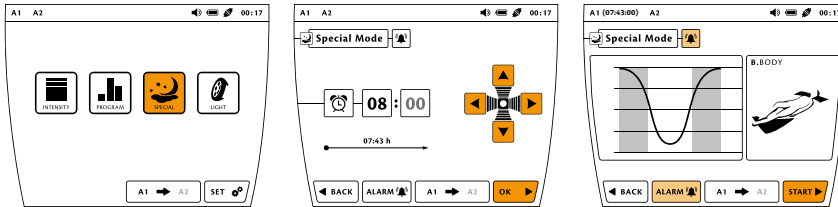
Special application – Sleep Program

- >> Press the button [SPECIAL].
- >> Using the arrow keys, set the desired time for the program's end (wake up time).
- >> Now activate the alarm signal if desired, by pressing the [ALARM] button.

- >> Confirm your selection by pressing the [OK] button.
- >> Now start the Sleep Program with [START].

Please take into consideration that the Sleep Program requires a sleeping period of at least 6 hours and not longer than 10 hours, as it supports the regenerative processes that occur during the core sleeping phase. 

Example – Sleep Program (2 application modules connected, the program runs via port A1 using the application module B.BODY)

Light therapy

- >> Press the [LIGHT] button.
- >> Select the desired application duration by pressing the corresponding button [L1] – [L4].
- >> Select the desired light intensity by pressing the corresponding button [1-10].
- >> Now start the light therapy with [START].

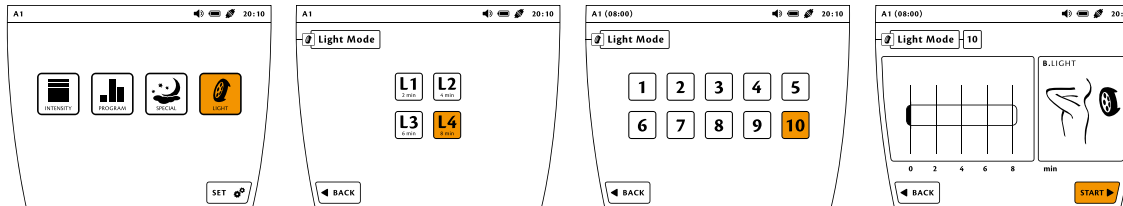
Usually high intensity levels are recommended. However for those with very sensitive skin, lower intensity levels can also be selected.

Do not let the LED light shine directly in the eyes. Use the provided protective goggles!

Example – light therapy (application module B.LIGHT connected on port A1, duration: 8 min.)

Creating a program

- >> Press the [SET] button.
- >> Then press the [EDITOR] button.
- >> Now choose the desired storage location for your new program (buttons 4-12).

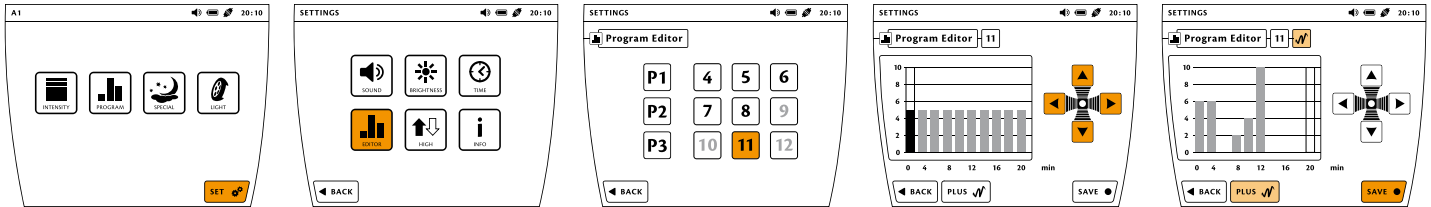
Procedure for creating programs

There are 10 blocks available. Each block has a duration of 2 minutes. Set the desired intensity by using the up or down arrow keys, select one of the blocks by using the left or right arrow keys.

>> If your program uses the *plus* Signal, press the [PLUS] button.

>> Now save the program with [SAVE].

Example – creating your own program

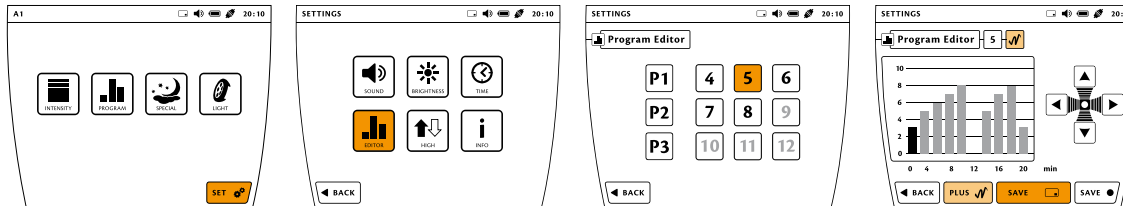
Saving on your memory card

- >> Slide the memory card into the card slot.
- >> Press the [SET] button.
- >> Press the [EDITOR] button.
- >> Now select the storage location of the desired program (using buttons 4-12).

- >> Make the necessary changes on the chosen program (see “Creating a program”).
- >> Now save the program on the memory card by pressing the button [SAVE – with the card symbol].

Example – saving on a memory card

Default settings – application parameters

The parameters [PLUS] and [INTENSITY] can be adjusted to be the default setting.

Setting the *plus* Signal as default

- >> To switch the *plus* Signals on / off, press the [PLUS] button before starting the program.

Default setting signal intensity (only for B.PAD)

While traveling, the B.PAD can be used as a substitute for the B.BODY by reducing the intensity to [LOW].

- >> Press the [SET] button.
- >> Press the button [HIGH / LOW] button.
When pressed, the button switches from [HIGH] to [LOW] and vice versa.

Basic settings – B.BOX Professional

Setting the signal volume

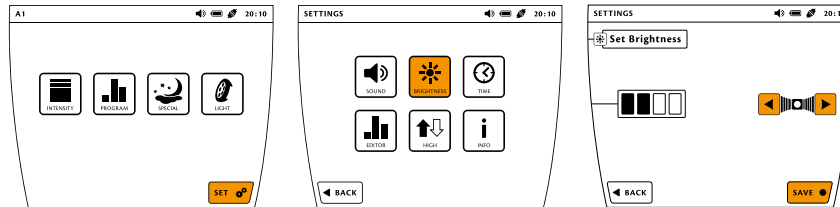
- >> Press the [SET] button.
- >> Press the [VOLUME] button.
- >> Set the desired signal volume with the help of the arrow keys.
The selection is always acknowledged by a signal in the chosen volume.
- >> Press the [SAVE] button to save your settings.

Confirmation of the [SET] button is possible only if no application is active.

Adjusting the display brightness

- >> Press the [SET] button.
- >> Press the [BRIGHTNESS] button.
- >> Now set the brightness with the help of the arrow keys. The display appears in the desired brightness.
- >> Press the [SAVE] button to save your settings.

Example – setting the display brightness

Setting the time

- >> Press the [SET] button.
- >> Press the [TIME] button.
- >> Now set the current time with the help of the arrow keys.
- >> Now press [SAVE] to save your settings.

To access the system information

- >> Press the [SET] button.
- >> Press the [INFO] button.
- >> Press the [OK] button to close the [INFO] window.

Updates

BEMER devices are a secure investment and are therefore designed to transfer product advancements with the corresponding updates to the B.BOX through the communication or maintenance interface. Since the nature, size and scale make different approaches for this possible, you will receive exact instructions for proceeding with any updates along with the delivery.

MAGYAR

FRANÇAIS

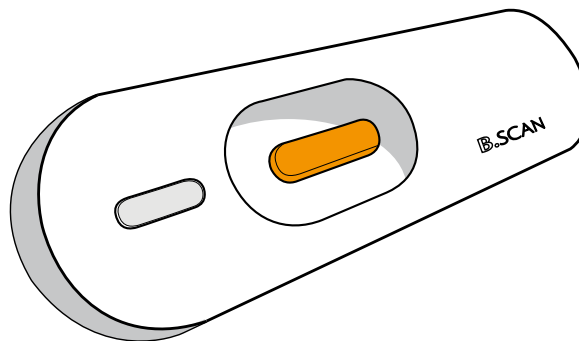
ENGLISH

DEUTSCH

B.SCAN

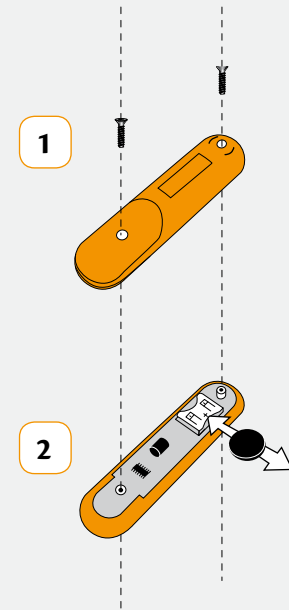
The B.SCAN is integrated in your B.BOX and is designed to check your BEMER devices along with safety precautions for their proper function. The B.SCAN makes the multidimensional BEMER Signal configuration audible and is therefore also ideally suitable as an effective display demonstration tool for interested persons.

The device is switched on by pressing the button once; pressing and holding down the button (for several seconds) will switch off the device. If no Signal system is detected, the B.SCAN automatically switches off after 10 seconds.

Changing the battery of the B.SCAN

- >> Remove both screws (see Fig. 1).
- >> Now carefully remove the back panel.
- >> Now slide the old battery sideways (in the direction of the arrow) away from its holder (see Fig. 2).
- >> Now slide the new battery (type CR 2032) into the holder. Please note the “+ pole” battery symbol (see Fig. 2).
- >> Now place the back panel back on the device and screw it back together.

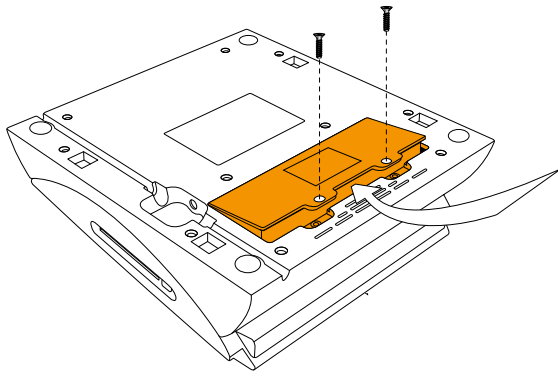
Rechargeable battery – B.BOX (optional accessory)

The high capacity, rechargeable battery used with the B.BOX makes it independent of electricity obtained from a power outlet and allows mobile application of the BEMER technology.

- >> Completely charge the battery before using for the first time.
- >> It takes about 6 hours to completely charge the unloaded battery.
- >> A newer battery, or one that has not been used over a period of time, will not be able to initially attain its full loading capacity.

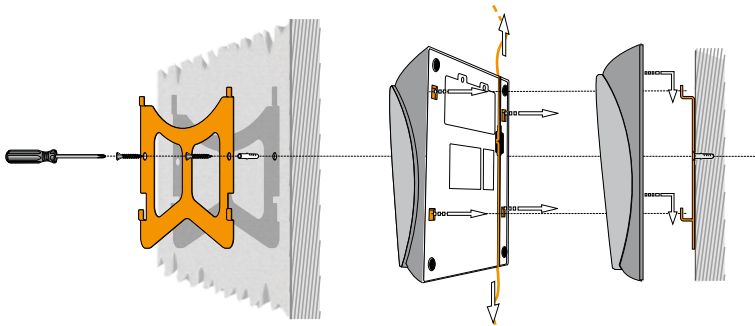
Assembly / Changing the battery

- >> Loosen the screws and remove the cover of the battery compartment.
- >> Tilt the rechargeable battery into the battery compartment and make sure that it is securely in place.
- >> Now attach the battery.

Wall mount assembly – Pro-Set (Classic-Set: optional accessory)

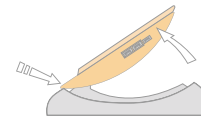
Assemble the wall mount as shown. Please note the course of the cable duct and use to assemble according to the circumstances.

Basic position recommended when using the wall mount.

Position 1 and position 2, recommended for table application with different angles.

MAGYAR

FRANÇAIS

ENGLISH

DEUTSCH

Cleaning

- >> Disconnect the power supply from the mains power supply.
- >> Clean the BEMER devices and its components with a dampened, non-abrasive cloth. Use only mild, non-abrasive, commercially available cleaning detergent (dishwashing liquid), according to the manufacturers instructions.
- >> All BEMER devices and all related components are neither waterproof nor rain-proof! **Do not allow moisture to penetrate into the device!**
- >> Do not use any abrasive cleaners, alkaline solutions or alcohol-based cleaners!
- >> Do not use any sharp edged or pointed objects.

Maintenance

BEMER products, when used correctly, are maintenance-free. Nevertheless, check your device condition regularly.

- >> Repair work may only be carried out by the manufacturer or by an authorized specialist, otherwise the guarantee is no longer valid.
- >> Before you send a BEMER product for repair to the manufacturer or to one of their authorized specialists, make sure to first carry out a generalized error troubleshooting check.
- >> For simplified processing of repairs or to rule out an error, contact the corresponding service address.
Use the original packaging for safer delivery.

Proper disposal

Properly dispose of the device or batteries according to local laws. You can receive further information from your sales dealer or directly from the manufacturer.

Dispose of used batteries correctly, do not throw in an open fire.

Danger of explosion!

Batteries do not belong in household garbage!

Troubleshooting

General information

BEMER devices are equipped with intelligent electronics. Therefore they can automatically detect and display errors. Should an error be displayed or should you suspect malfunction, then check the following:

- >> Is the electricity connected to the power supply (does the power supply have a corresponding control light)? If not, then connect the electricity to the power supply.
- >> Is the power supply connected to the B.BOX? If not, then connect the electricity to the power supply.
- >> Are the desired application modules properly connected?
- >> Are all the cables and plugs that come into question undamaged?

If an error still persists or if one has not been properly located, please follow the instructions of the corresponding error code or contact your responsible service address. When in doubt, you can get information at **WWW.BEMERGROUP.COM**, where you will find your next closest contact person.

Switching off the error message

When the error or cause has been remedied, you can switch off the error message on the device again. This procedure using the B.BOX Pro is self explanatory and is found in the menu navigation. The error message on the B.BOX Classic is shut off by pressing the button "Set".

ERROR MESSAGE	CAUSE	TROUBLESHOOTING
ERROR 1	An active applicator has been removed during the application.	Reconnect the applicator and start the application again.
ERROR 2	The temperature of the device is too high.	Check the environmental conditions and let the device cool down. When in doubt, contact service.
ERROR 3-4	There is a problem with the application module.	Please contact service.
ERROR 5-6	There is a malfunction of the power supply.	Please check for the possible sources of error such as: the power supply / the vehicle connection cable (12 volt) or the rechargeable battery (see technical data). When in doubt, contact service.
ERROR 7-8	There is a malfunction with the control device.	Please contact service.
ERROR 9-12	There is a malfunction of the application module.	Please contact service.
ERROR 101	There is no applicator connected.	Please connect an applicator to the device.
ERROR 104	An unknown applicator has been detected.	The applicator module is invalid, please contact service.
ERROR 105	The capacity of the battery is too low.	The battery is empty, recharge the battery by connecting to the power supply.

Technical Data

B.BOX Classic	Art. No. 420100
Dimensions (W x H x D)	22 x 23 x 7 cm
Weight	1.3 kg
Operating voltage	12 – 15 V DC
Power input	2 A
Performance	max. 30 W
>> Connections for 2 application modules	
>> Maintenance interface	

B.BOX Pro	Art. No. 420200
Dimensions (L x W x H)	22 x 23 x 7 cm
Weight	1.4 kg
Operating voltage	12 – 15 V DC
Power input	2 A
Performance	max. 30 W
>> Connections for 2 application modules	
>> Maintenance interface	
>> Communication interface	

B.BODY Classic	Art. No. 430100
Dimensions (L x W x H)	180 x 60 x 2 cm
Weight	2 kg
Surface material	PES
Coils	6
Average flux density	≈ 35 μT (max. level)
Average flux density <i>plus</i>	≈ 50 μT (max. level)

B.BODY Professional	Art. No. 430200
Dimensions (W x H x D)	180 x 60 x 2 cm
Weight	2 kg
Surface material	Alcantara, PES
Coils	6
Average Flux density	≈ 35 μT (max. level)
Average flux density <i>plus</i>	≈ 50 μT (max. level)

Vehicle connection cable	Art. No. 440300
Primary voltage	12 V

Rechargeable battery	Art. No. 440200
Dimensions (W x H x D)	14 x 6 x 2 cm
Weight	220 g
Operating voltage	7.2 V
Battery type	Li-Ion
Capacity	4500 mAh
Operating time	50 applications *

B.BOX power supply	Art. No. 440100
Primary voltage	100 – 240 V AC
Primary current	1 A
Frequency	50 – 60 Hz
Secondary voltage	15 V DC
Secondary current	2.5 A max.
Protection class	SK II
Primary connection	1 m
Secondary connection	1.8 m

>> Power supply	EU cable, US cable, UK cable, AU cable
Power supply indicator	LED green

B.SCAN	Art. No. 450100
Dimensions	12 x 3 x 2 cm
Weight (incl. battery)	39 g
Power supply	Battery CR2032
Operating time	≈ 1000 Measurements
Power input	Approx. 2 mA

B.GRIP	Art. No. 431000
Dimensions (W x H x D)	12 x 12 x 5 cm
Weight	200 g
Connection cable	2 m
Connection plug	6 pin

* Value are taken from the measurements (by a respective short application with the level 5) with a connected application module.

B.SPOT	Art. No. 431100
Dimensions (W x H x D)	13 x 13 x 3 cm
Weight	300 g
Coils	1
Average flux density	≈100 μT (max. level)
Average flux density <i>plus</i>	≈100 μT (max. level)

B.LIGHT	Art.No. 431200
Dimensions (W x H x D)	13 x 13 x 3 cm
Weight	200 g
Light emitting diodes (LED)	52 pieces SMD with spec. devel. lens
Wavelengths	660 nm
Intensity	4000 mcd

B.BODY Attachment belt	Art. No. 450800
Dimensions (L x W x H)	180 x 70 x 0.2 cm

Attachment belt (B.GRIP)	Art. No. 450700
Dimensions (L x W x H)	730 x 120 x 4 mm
Weight	50 g
Surface material	Neoprene

B.PAD	Art. No. 430300
Dimensions (L x W x H)	1110 x 130 x 15 mm
Weight	1000 g
Coils	3
Average flux density	≈100 μT (max. level)
Average flux density <i>plus</i>	≈150 μT (max. level)
Surface material	Neoprene

Protective goggles	Art. No. 450600
Glazing:	Shade 2, GA 166 CE (anti-fog, scratchproof, 100% UV protection), impact resistant polycarbonate. Do not use goggles 5 years after the date of purchase.

Environmental conditions during operation

Temperature	+10...+40 °C
Air pressure	700...1060 hPa \cong 4000 m (ASL)
Relative humidity	30...75% RH, non-condensing

Storage and transportation

Temperature	-20...70 °C
Air pressure	500...1060 hPa \cong 4000 m (ASL)
Relative humidity	10...75% RH, non-condensing

Standard marking

EN 60601-1	Medical electronic equipment
EN 60601-1-2	Electromagnetic compatibility
EN 60601-1-4	Programmable electronic medical systems
EN 60601-1-6	Usability
EN ISO 14971	Application of risk management on medical products
EN 62304	Software life cycle processes
EN 62366	Application of usability on medical products

Guideline and manufacturer's declaration of conformity (tables according to EN 60601-1-2)

Electromagnetic emission

The BEMER SET CLASSIC / PRO is designed for operation in the electromagnetic environment listed below.

The customer or user of the BEMER SETS CLASSIC / PRO should ensure that the device be used in such an environment.

EMISSION TESTS	COMPLIANCE	ELECTROMAGNETIC ENVIRONMENT – GUIDELINES
RF emissions according to CISPR 11	Group 1	The BEMER SET CLASSIC / PRO uses RF energy exclusively for its internal function. Therefore, its RF emission is very small. It is unlikely that nearby electronic devices will be disturbed. The BEMER SET CLASSIC / PRO is designed for use in all premises, including such residential areas that are directly connected to a public power supply that also supplies residential buildings.
RF emissions according to CISPR 11	Class B	
Harmonics according to IEC 61000-3-2	Not usable	
Voltage fluctuations / flicker according to IEC 61000-3-3	Not usable	

Table according to EN 60601-1-2

Recommended safety distance between portable and mobile RF telecommunication devices and the BEMER SET CLASSIC / PRO

The BEMER SET CLASSIC / PRO is designed for operation in an electromagnetic environment, in which the RF disturbances are controlled. The customer or the user of the BEMER SETS CLASSIC / PRO can help to prevent electromagnetic disturbances by keeping a minimum distance between portable and mobile RF telecommunication devices (transmitter) and the BEMER SET CLASSIC / PRO – depending on the power output of the communication device, as indicated below.

OUTPUT POWER OF THE TRANSMITTER W	SAFETY DISTANCE, DEPENDENT ON THE TRANSMITTED FREQUENCY m		
	150 KHZ TO 80 MHZ $d=1.2\sqrt{P}$	80 MHZ TO 800 MHZ $d=1.2\sqrt{P}$	800 MHZ TO 2.5 GHZ $d=2.3\sqrt{P}$
0.01	0.12	0.12	0.23
0.1	0.38	0.38	0.73
1	1.2	1.2	2.3
10	3.8	3.8	7.3
100	12	12	23

For transmitters whose maximum output power rating is not listed in the table above, the distance can be determined by calculating the equation on each column, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

Electromagnetic immunity

The BEMER SET CLASSIC / PRO is designed for operation in the electromagnetic environment listed below.

The customer or user of the BEMER SETS CLASSIC / PRO should ensure that the device is used in such an environment.

EMISSION TESTS	IEC 60601 TEST LEVEL	COMPLIANCE LEVEL
Electrostatic discharge according to IEC 61000-4-2	± 6 kV Contact discharge ± 8 kV Air discharge	± 6 kV Contact discharge ± 8 kV Air discharge
Fast transient electrical interferences according to IEC 61000-4-4	± 2 kV for power supply lines	± 2 kV for power supply lines
Surge voltage / surges according to IEC 61000-4-5	± 1 kV Voltage outer conductor – outer conductor ± 2 kV Voltage outerconductor – earth	± 1 kV Voltage outer conductor – outer conductor ± 2 kV Voltage outerconductor – earth
Voltage dips, interruptions and variations of the power supply according to IEC 61000-4-11	< 5 % U_T for 1/2 cycles (> 95 % dip in U_T) < 40 % U_T for 10 cycles (60 % dip in U_T) < 70 % U_T for 25 cycles (30 % dip in U_T) < 5 % U_T 5 s (> 95 % dip in U_T)	< 5 % U_T for 1/2 cycles (> 95 % dip in U_T) < 40 % U_T for 5 cycles (60 % dip in U_T) < 70 % U_T for 25 cycles (30 % dip in U_T) < 5 % U_T 5 s (> 95 % dip in U_T)
Power frequency (50 / 60 HZ) magnetic field according to IEC 61000-4-8	3 A/m	3 A/m

Note: U_T is the alternating power supply voltage prior to application of the test level

ELECTROMAGNETIC ENVIRONMENT – GUIDELINES

	The quality of the power supply should correspond to that of a typical commercial or hospital environment.
	The quality of the power supply should correspond to that of a typical commercial or hospital environment.
	The quality of the power supply should correspond to that of a typical commercial or hospital environment.

Electromagnetic immunity

Portable and mobile radio equipment should be used no closer to the BEMER SET CLASSIC / PRO and its power cables, than the recommended safety distance calculated from the equation applicable to the frequency transmitter.

EMISSION TESTS

IEC 60601 TEST LEVEL

COMPLIANCE LEVEL

Conducted RF interference according to IEC 61000-4-6	3 V _{actual value} 150 kHz to 80 MHz outside of the ISM bands	3 V	
Radiated RF interference according to IEC 61000-4-3	3 V/m 80 MHz to 2.5 GHz	3 V/m	

NOTE 1: At 80 MHz and 800 MHz, the higher the frequency range applies. NOTE 2: These guidelines may not apply in all situations. The propagation of electromagnetic waves is affected by the absorption and reflection from structures, objects and people.

a) Field strength from fixed transmitters such as base stations for radio (cellular / wireless) telephones and land mobile radios, amateur radio, AM / FM radio stations and TV stations can not be predicted theoretically with accuracy. To determine the electromagnetic environment of a stationary transmitter, a site survey should be considered. If the measured field strength in the location in which the BEMER SET CLASSIC / PRO is used exceeds the conformity level listed above then the BEMER SET CLASSIC / PRO should be monitored to verify normal operation. If abnormal performance is observed, additional measures might be necessary, such as readjusting the BEMER SETS CLASSIC / PRO or changing to another location.

ELECTROMAGNETIC ENVIRONMENT – GUIDELINES

Recommended safety distance:

$$d = 1.2\sqrt{P}$$

$$d = 1.2\sqrt{P} \text{ for 80 MHz to 800 MHz}$$

$$d = 2.3\sqrt{P} \text{ for 800 MHz to 2.5 GHz}$$

Where P is the output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended safety distance in meters (m). The field strength from a stationary radio transmitter, according to an onsite survey ^a is less than the compliance level ^a for all frequency ranges.

b) The field strength is less than 3V/m over the frequency range of 150 kHz to 80 MHz.

Guarantee

The manufacturer guarantees that the BEMER products, when used under normal operation and maintenance conditions, perform error-free and are made of trouble-free material quality for a period of 36 months from the date of delivery.

This guarantee does not extend to cover possibly required recalibration and / or maintenance of any kind.

The guarantee applies only to the first acquiring purchaser and does not extend to products or components thereof, which were not properly used or to products where changes have been made. Any guarantee is rendered null and void by non-intended use of the BEMER products.

The guarantee obligation is limited to repairs or the exchange of a product that has been sent to the manufacturer during the guarantee period. A prerequisite is that the manufacturer acknowledges that the product is defective and the defect is not attributable to improper handling or changes to the device or due to any abnormal operating conditions.

Any guarantee obligation expires, should repairs on BEMER products be carried out by persons not authorized from the manufacturer.

The terms of this guarantee apply exclusively and in place of all other contractual or statutory obligations, including, but not limited to the legal warranty regulations of the marketability, the suitability and the expediency for the specific application.

The manufacturer assumes no responsibility for direct and indirect, incidental or consequential damages, regardless of whether they are the results of lawful, unlawful or any other actions.

Guarantee and liability claims for personal and property damage are excluded, should they result from one or more of the following reasons:

- >> Unintentional or negligent use of BEMER products.
- >> Improper use, operation and maintenance of BEMER products.
- >> Noncompliance with the guidelines / warnings found in the operating instructions.
- >> Unauthorized, structural changes on BEMER products.
- >> Noncompliance with the operating instructions.
- >> Improperly handled repairs and the use of parts, replacement parts or accessories not recommended by the manufacturer.
- >> Disaster and force majeure.

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BEMER Int. AG